



## Sartolab® Multistation

Hands-Free Parallel Filtration of up to Six Samples

### User Benefits

- One single vacuum source to enable the filtration of up to 6 samples at the same time
- Hands-free filtration
- Time saving  
(no extra installation time to vacuum for each single filter unit)



### Product Information


Sartolab® MultiStation is a stand especially designed to hold 1 to 6 vacuum filtration units. It allows a filtration of all samples at the same time. The MultiStation is permanently connected to a vacuum source; the vacuum filtration units can be easily installed in the MultiStation for quick and easy filtration of samples without extra time consuming connection.

Sartolab® MultiStation functions with the vacuum filtration units Sartolab® RF 50. The adapter connecting the funnel to the 50 mL sample tube perfectly suits into the bracket of the MultiStation.

One click and the filtration unit hangs in the device. That assures a perfect stability of the filter; with a second click, the vacuum is automatically applied and the filtration begins.

The rotating head is multi-directional and enables easy handling of your samples in front of you. Furthermore, all 6 brackets are numbered to avoid any confusion during the filtration of multiple samples.

## Example of Filtration of 6 x 50 mL Samples with and without the Sartolab® MultiStation

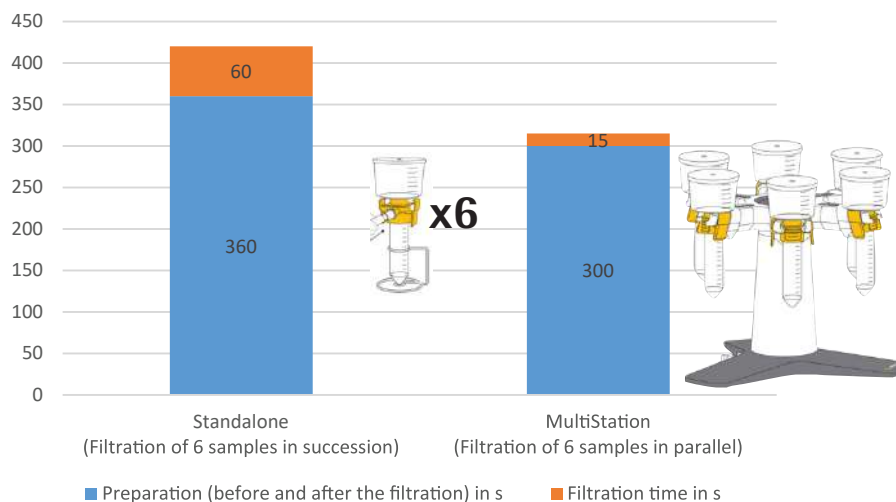

## Specifications

Material (visible parts)	ABS/Aluminum/Stainless steel
Dimensions (Length x Width x Height)	307 x 348 x 281 mm
Weight	4.6 kg
Tubing adaptor	Designed for tubing with an inner diameter between 4 and 10 mm and with a wall thickness of minimum 3 mm.

## Ordering Information

Order Number	Description	Quantity
SDLC01	Sartolab® MultiStation for the filtration of up to 6 samples in parallel	1 unit

## Consumables

Order Number	Description	Quantity
180E01-----2	Sartolab® RF 50, 0.22 µm, PES, vacuum filtration units for 50 mL samples	24 units
180F01-----2	Sartolab® RF 50, 0.45 µm, PES, vacuum filtration units for 50 mL samples	24 units

**Sartorius Lab Instruments**  
GmbH & Co. KG  
Otto-Brenner-Strasse 20  
37079 Goettingen, Germany

Phone +49.551.308.0  
Fax +49.551.308.3289  
www.sartorius.com

USA Toll-free +1.800.635.2906  
UK +44.1372.737159  
France +33.1.70.62.50.00  
Italy +39.0362.5557.11  
Spain +34.913.586.095  
Russian Federation +7.812.327.53.27  
Japan +81.3.3740.5408

Specifications subject to change without notice.  
Copyright Sartorius Lab Instruments GmbH & Co. KG.